DENTAL UNIT CESTAI Rod Type OPERATING INSTRUCTIONS

IMPORTANT

- This manual provides operating instruction for CLESTA II ROD TYPE.
- The instructions contained in this booklet should be thoroughly read and understood before operating the unit.
- After the installation is completed, file this manual and refer back to it for future maintenance.
- If you have any questions about this Manual or this product, please contact us.
 If manual becomes unreadable or is lost, please request a new manual by contacting your dealer.
- Installation should be conducted by authorized personnel only. Follow instructions on installation manual.



TABLE OF CONTENTS

	Page
SAFETY PRECAUTIONS	1
1. OVERVIEW AND MAJOR PARTS	6
2. DIMENSIONS & SPECIFICATIONS	
2-1. DIMENSIONS	7
2-2. SPECIFICATIONS	7
3. OPERATING INSTRUCTIONS	
3-1. DOCTOR TABLE SECTION	8
3-2. CUSPIDOR UNIT SECTION	15
3-3. FOOT CONTROL SECTION	18
4. SAFETY LOCK DEVICE	18
5. CARE AND MAINTENANCE	19
6. MAINTANANCE AND INSPECTION	
6-1. GUIDELINE FOR DAILY MAINTENANCE AND INSPECTION	25
(MAINTENANCE AND INSPECTION BY USER)	
6-2. GUIDELINE FOR PERIODICAL CHECK-UP	27
7. DETACHABLE PARTS	27
8. ELECTROMAGNETIC COMPATIBILITY	28
9. HANDPIECES THAT CAN BE USED WITH THIS PRODUCT	30
10. CHAIRS THAT CAN BE USED WITH THIS PRODUCT	30
11. DENTAL LIGHTS THAT CAN BE USED WITH THIS PRODUCT	30
12. ACCESSORIES	31

Intended Use of the Product

This product is an active therapeutic device intended for the exclusive use for diagnoses, treatments and relative procedures of dentistry.

The product must be operated or handled by the qualified dentists or by dental staffs under the supervision of the dentist.

Such dentists or dental staffs should instruct and/or assist the patients to approach to and leave from the product.

Patients should not be allowed to operate or handle the product unless he/she is so instructed. The product is supplied together with the handpieces like electric micromotor, air turbine and/or motor, scaler and so on.

In addition, it is the device which combines with dental chair and dental light.

Compliance with Directives

CLESTAII Compliance with the following directives, Medical Device Directive 93/42/EEC, RoHS Directive 2011/65/EU.

Handpieces that can be used with this product

Use the Handpieces described on page 30.

Chair that can be used with this product

Use the Chair described on page 30.

Dental Light that can be used with this product

Use the Dental Light described on page 30.

Important Notes

In case of the troubles, please contact Takara Belmont offices or your dealers.

Do not disassemble or attempt to repair.

Disassembly, repair or modifications shoud only be done by a qualified repair technician.

Attempts at disassembly, repair or modifications may lead to abnormal operation and accidents.

Disposal of residue material

Please request a special contractor when you dispose amalgam.

In case of disposal of equipment

When disposing the chair, appropriately dispose complying with all current applicable regulations and local codes.

In EU area, EU directive 2002/96/EC on waste electrical and electronic equipment (WEEE) is applied on this product. In this directive, environment conscious recycling/abandonment is obligated.

SYMBOLS

In this manual, on the labels or on the control panel of CLESTA II ROD TYPE, following symbols are used. Confirm the meaning of each symbol.

	Protective earth (ground)	Ļ	Functional earth (ground)		ON (power)	\bigcirc	OFF (power)
LP	Chair last position	0	Chair auto return	1	Chair preset 1	2	Chair preset 2
1	Chair auto control	Ś	Chair manual control	\wedge	To raise the chair		To recline the backrest
\checkmark	To lower the chair	Z	To raise the backrest	-Ŏ.	Fiber optic handpiece light on//off	Ĭ	Handpiece coolant spray on/off
	Rotation mode select	Ē	Micro motor Forward/Reverse select	F	Function	€	Store
Å	Rotation speed contol	$\mathbb{N}^{\mathbb{A}}$	Scaler power control	$ \leq$	Syringe	ζ''	Bowl flush
Щ	Cupfiller	<u> </u>	Dental light on/off		Dental light mode selection	min.	Minus
sec.	Plus	\ / 7 F	Service outlet (water)	עלאיד	Service outlet water flow control	\\/ ᠯᠮ	Service outlet (air)
W	Water	A	Air		Water heater	\sim	alternating current
	General warning sign		Caution	×	Type B Applied Parts	REF	Reference number
IPX 1	Classification of foot controller		Refer to operating instructions		Manufacturer	~	Date of manufacture
SN	Serial number	CE 0197	Comply with MDD,93/42/EEC	X	Separate collection for electrical and electronic equipment	EC REP	Authorized representative in the European community
135℃ ∫∫∫∫	Autoclave Symbol This symbol on component means that the component can be sterilised with an autoclave at 135°C max.						

*1 : The color of black part is blue in the actual label.

REQUEST FOR WATER FLUSHOUT (DISCHARGE)

Precautions for water quality safety

E type:Standard A type:Option Ε

А

To maintain the water quality for treatment and for stable water supply to the handpieces, ensure to carry out flushout of the residual water in the unit before starting treatment.

At the beginning of the week or in the morning, that is, a long period of idle time, the residual water in tubings and water heater in the unit tends to permit bacteria being propageted.

Ensure to carry out flushout of water lines in the unit before starting treatment, for the safe treatment and for the operation of the handpiece without any problem.

It is also recommended to flush out the residual water in the unit and in handpiece tubings with fresh water, after the treatment in the morning and also in the evening to prevent bacteria propagation.

Standard length of time required for flushout of water lines in the unit			
Handpiece line	About 40 seconds per each of turbine, motor, scaler and syringe.		
Air Turbine	(About 40 seconds also in case of flushout of all of them in a batch.)		
Micro Motor	ACAUTION		
Scaler	Apply flushout to the scaler with the body kept fitted. Failure		
Syringe (for both doctor / assistant)	may result otherwise.		
Cuspidor line			
Cupfiller (water heater)	About 5 minutes for cleaning water in cupfiller line.		
Bowl flush			

Flushout procedure E		
Handpiece line Pick up one handpiece out of the holder, and hold it over the cuspidor bowl with the turbine kept as it is and the contra angle connected to the motor. (In case the flushout function is provided, pick up all the handpieces in a lump from the holder, and hold them on the cuspidor bowl.)	Short time flush out: For Flushing out the Handpiece water line only As for flushout procedures, please refer to Page 12 selection and operation of short time flush out.	Long time flush out: For Flushing out the cupfiller and cuspidor bowl after handpiece flushout As for flushout procedures, please refer to Page 12
Cuspidor line Cupfiller (water heater) Bowl flush		Selection and operation of long time flush out.

Flushout procedure A	
Handpiece line Pick up one handpiece out of the holder, and hold it over the cuspidor bowl with the turbine kept as it is and the contra angle connected to the motor. (*In case the flushout function is provided, pick up all the handpieces in a lump from the holder, and hold them over the cuspidor bowl.)	* Applicable when flush-out function was selected as optional. Turn the toggle switch on and flush-out the handpiece
Cuspidor line Cupfiller (water heater) Bowl flush	Repeat cupfilling a few times. (at least 3 times) Press the bowl flush switch and keep flushing the cuspidor bowl for 5 minutes.

■ Before use, read the "Safety precautions" carefully to ensure proper use.

The following information is designed to ensure safe use of this product and to prevent injury and damage to you and others. The precautions contained here are classified depending on the severity and degree of imminence of possible injury or damage resulting from improper use. Be sure to follow all the information, which is important for safety.

Classification of precautions	Severity and degree of imminence of possible injury or damage
WARNING	This symbol indicates that "ignorance of these precautions may lead to severe injury or even death as a result of improper use."
CAUTION	This symbol indicates that "ignorance of these precautions may lead to mild or moderate physical injury or damage to property as a result of improper use."
NOTICE	This symbol indicates that "it is recommended to follow these precautions for safety."

WARNING

1. Be sure to turn off breakers for equipment in the clinic when this product will not be used for a long period of time

Be sure to turn off breakers for equipment in the clinic when this product will not be used for a long period of time (following the completion of work, during the suspension of business, etc.). Insulation degradation may cause electrical fire.

2. Be sure to turn off the main switch upon completion of work or during work breaks

Be sure to turn off the main switch upon completion of work or during work breaks. This prevents incorrect operation due to accidental contact and associated hazards.

3. Do not sit on other than seat

When the backrest is at the forward position. do not sit on or place an undue load on the headrest or legreat of dental chair. This could cause the unit to topple or could damage the unit.

4. Do not place an undue load on the arm

Do not get on or place an undue load on the arm of this unit or dental chair armrest. This could cause the unit to topple or other accidents.

5. Be sure to establish a grounding connection

Be sure to establish a proper grounding connection. (Refer to a vendor for grounding connection.) Failure or electric leakage may lead to electric shock.

6. Never disassemble, repair or modify this product

Individuals other than certified repair technicians should not disassemble or attempt to repair and modify this product. This could lead to an accident, failure, electric shock or fire.

7. Use with caution in the presence of electromagnetic interference waves

Do not place this product around equipment generating electromagnetic waves (including communications equipment, elevators, etc.) as incorrect operation of this product may occur in the presence of electromagnetic interference waves. Do not use equipment generating electromagnetic waves, such as mobile phones, around this product.

8. Be sure to turn off the main switch when electrocautery is in use

Be sure to turn off the main switch when electrocautery is in use, because noise may cause incorrect operation of this product.

WARNING

9. Ensure the maintenance of this product

• Failure to maintain this product may lead to physical injury or property damage.

• Refer to the section of maintenance.

10. Do not place objects weighing 1.5 kg or more on the Sub tray

Do not place objects weighing 1.5 kilograms or more on the Sub tray. This could cause damage to the Doctor's table, defective function or accidents.

11. Be sure to use the mirror cover

• Be sure to use the mirror cover of the dental light when the light is turned on. Direct contact with lamps may cause burns.

• See the Instruction Manual of the dental light for further information.

12. Be sure to turn off the power when replacing lamps

• Be sure to turn off the power when replacing the dental light. This could result in electric shock.

• Use only dedicated halogen lamps.

• Immediately after a halogen lamp has burnt out, the lamp and the lamp holder are still hot. Replace the lamp after they cool down.

• Do not touch halogen lamps with bare hands.

• See the Instruction Manual of the dental light for further information.

13. Immediately wipe off any water spills or leakage on the floor

Immediately wipe off any water spills or leakage on the floor. Decreased strength of the floor may lead to physical injury including fall, or property damage.

14. Use with caution on patients with a cardiac pacemaker

Use this product with extreme caution on patients with a cardiac pacemaker. In the case of any abnormalities in patients during use, immediately turn off this product and discontinue use.

15. Prohibition of maintenance

During operation, repair and maintenance are prohibited.

16. Prohibition of using this equipment adjacent to or stacked with other electronic equipment

Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

17. Prohibition of placing portable RF communications equipment adjacent to this product

Portable RF communications equipment (including peripherals such as antenna cablesand external antennas) should be used no closer than 30 cm (12 inches) to any part of the CLESTAII, including cables specified by the manufacturer. Otherwise, degradation of the performance of this euipment could result.

ACAUTION

1. Only experienced personnel should use this product

Only dentists or other dental professionals should use this product.

2. Confirm safety before use.

Before use, confirm that the parts are correctly and safely operating and that there are no obstacles around this product.

3. Pay attention to patients and children

Keep your eyes on patients (especially, children) so that mischief or inadvertent operation of equipment will not lead to unexpected accidents.

A CAUTION

4. Discontinue use if you feel that "something is wrong"

Always be careful to inspect this product for looseness, rattling, tilting, wobbling, sounds, temperature, odors, etc. Immediately discontinue use at the first feeling that "something is wrong."

5. keep your eyes on the patient during operation.

• Confirm that the patient is seated in the proper position. Keep your eyes on the patient during the operation.

• Pay special attention to surroundings at automatic operation of the dental treatment table. Damage to the backrest, stool or Doctor's table may occur.

6. Pay attention during movement of the Doctor's table

• Pay attention to surroundings when you move the Doctor's table. Injury by the tips of handpieces, etc., may occur.

• Be sure to move the Doctor's table by holding the handle of the unit.

• Be sure to move up & down the Doctor's table by releasing the balance arm brake of the unit.

7. Do not place anything hot on the unit

Do not place anything hot on the unit. This could cause deformation or discoloration.

8. Do not smack or rub this product

Do not smack or rub this product forcefully. This could cause damage to covers or defective function.

9. Precautions for cleaning a spittoon bowl

The spittoon bowl is made of glass. Handle with care. Do not wash it with hot water. Otherwise, it may break.

10. Pay attention to water discharge when attach/detach the spittoon bowl

Make sure to turn the main switch off when you attach/detach the spittoon bowl. Cup filler sensor may react to hands. Resulting water spill.

11. Precautions for cleaning the resin cover

For cleaning, do not use cleaning agents containing solvent or abrasives, thinners or oil-based alcohol (butanol and isopropyl alcohol), which may cause cracks.

12. Immediately wipe off drug solution when it comes into contact with this unit

Should drug solution or water comes into contact with this unit, immediately wipe it off with a dry soft towel, etc. This could result in defective function or electric leakage as well as spotting or rusting.

13. Close the water main valve and turn off the main switch upon completion of work

Be sure to close the water main valve and the main switch at the end of each work day to prevent water leakage from occurring.

14. Be sure to operate switches with your hands

Be sure to operate switches with your hands, except the foot controller, which is operated with your foot. Operation with body parts other than hands may cause damage or incorrect operation.

15. Pay attention during the headrest operation

Do not allow hands, fingers, or hair to become entangled in the moving parts of the headrest during operationg

16. Be careful not to heat an empty water heater

Exercise caution as heating of an empty water heater may result in burning of the heater, leading to fire.

17. Precautions when using handpieces, etc.

In order to ensure safety, be sure to confirm that rotation has completely stopped before changing the rotation speed of the micromotor or inserting/removing the bar of handpieces. See the Instruction Manual accompanying various handpieces for further information.

ACAUTION

18. Precautions for cleaning the operation panel (membrane switches)

Penetration of droplets of sanitizing spray into the back of the operation panel may be associated with switch failure. Use a paper towel soaked with sanitizing solution to clean the surface of the operation panel.

19. Precautions for cleaning

• Never use sandpaper, metal scrub brushes and abrasive cleaning agents to clean the unit.

• Do not use acidic cleaning agents or alkaline pipe cleaning agents to avoid corrosion of metal and resin parts, etc. Never use cleaning agents that include sodium hypochlorite.

20. Practice of flush out

Practice the flush out of water retained in the unit before the start of clinical practice at the beginning of each work day to maintain the quality of water for dental treatment and ensure a steady supply of water to handpieces. Refer to the procedure for flush out.

21. Set the pressure of the water tank at 200 kPa or less

Adjust the air supply pressure for the water tank to 200 kPa or less. An excessively high pressure may cause damage to the water tank.

22. Do not use water other than purified water, distilled water or pure water for the water tank

The water tank is intended only for use with purified water, distilled water and pure water. Do not use mouthwash or electrolyzed water, such as ConCool or povidone iodine, as they may cause clogged tubing or affect internal valves and equipment.

23. Read the documents accompanying the various pieces of equipment

Before use, be sure to carefully read the package inserts and Instruction Manuals accompanying the various pieces of equipment (including optional articles) to ensure proper use.

24. Precautions for sterilization

Do not sterilize other than specified method. It may cause device failure.

25. Combination with accessories or other devices

Use only our qualified accessories, other devices or equipment for this product to ensure the safety of product.

26. Pay attention to allergic reaction of a patient

While a vaccume tip is placed in a patient mouse, pay attention to allergic reaction. If allergic reaction start to happen, immediately stop the usage of a vaccume tip.

27. Precautions for sterilization

Do not sterilize other than specified method. It may cause device failure.

NOTICE

1. Troubleshooting and contact information

In the case of any problems, discontinue use, turn off the main switch and contact the dealer or our company.

2. Check operation of the compressor

With no air supplied, this product does not operate even after turning on the main switch. Turn on the power of the compressor before operating this product.

3. Precautions when using water other than tap water

The water unit is intended for use with tap water. Caution should be exercised as the use of water other than tap water (water through a sterilizer of water systems, etc.) may result in failure of equipment.

4. Use the turbine with a water check valve

Use the turbine with a water check valve. Contact the dealer or our company when a turbine without a water check valve will be used.

5. Handling of equipment in the case of a power failure

Put the handpiece in the holder and turn off the main switch if equipment stops working during use due to a power failure or other reasons.

1. OVERVIEW AND MAJOR PARTS

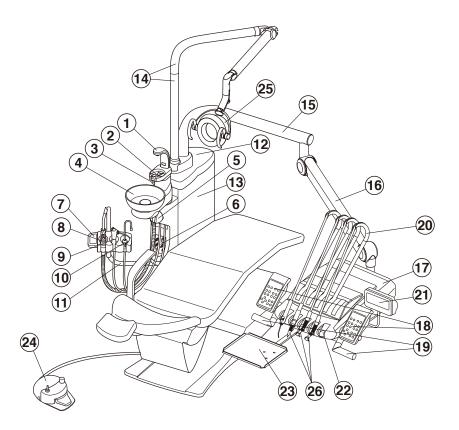


Fig.1-1 Overview and Major Parts

- 1 Cupfiller Nozzle
- 2 Cupfiller Base
- (3) Bowl Flush Nozzle
- (4) Spittoon Bowl
- (5) Solid Collector
- (6) Cuspidor Unit Control Panel
- Assistant Instruments Holder(With Membrane Switch)
- **8** Vacuum Handpiece
- (9) Assistant Side Syringe

- **10** Saliva Ejector Handpiece
- (1) Assistant Holder Arm
- (12) Cuspidor Top
- (13) Cuspidor Body
- (14) Light Pole
- (15) First Arm (Horizontal Arm)
- **16** Balance Arm
- (17) Table Housing
- (18) Main Control Panel
 - Right Hand Side : Standard Left Hand Side : Option

(19) Handle


Right Hand Side : Standard

- Left Hand Side : Option
- **20** Rod
- (21) Dental Size Film Viewer
- (22) Handpiece Rest
- **23** Sub Tray Holder
- **24** Foot Controller
- **25** Dental Light
- (26) Handpieces (Micromotor, Air Turbine/Motor, Scaler and etc.) Manufacturers recommend

to use the handpieces with CE markings

2. DIMENSIONS & SPECIFICATIONS

2-1. DIMENSIONS

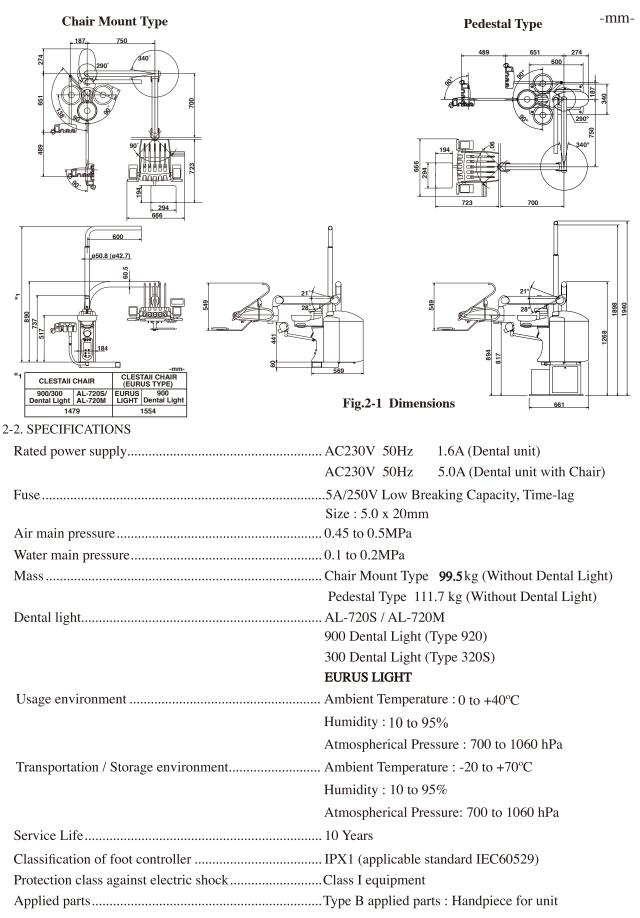

Equipment that is not suitable for use in air, flammable anesthetic gas, oxygen or nitrous oxide

3. OPERATING INSTRUCTIONS

3-1. DOCTOR TABLE SECTION

```
(1) Master Switch (Fig.3-1)
```

Turn on the master switch located on the right hand side (facing) under the doctor table, the power indicator on the main control panel illuminates in green.

Note : E Type

Fig.3-2 is indicated and it turns to Fig.3-3 after seveal seconds on the indicator. Wait operation until Fig.3-2 is changed to Fig.3-3 on the indicator.

(2) Main Control Panel



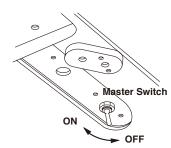


Fig.3-1 Master Switch and Indicator

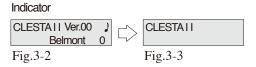

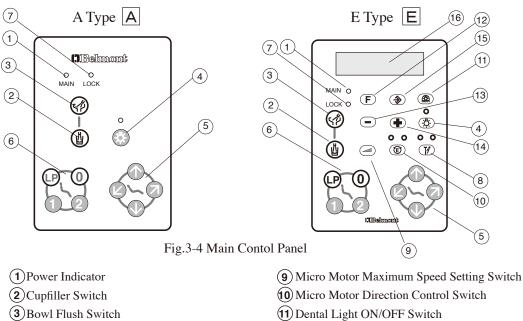

- (4) Light Pack Switch
- (5) Chair Manual Control Switch
- (6) Chair Auto Mode Control Switch
- (7)Safety Lock Indicator
- (8) Spray Water ON/OFF Switch

- (12) Function Switch
- (13) Decreasing Switch
- (14) Increasing Switch
- (15) Store Switch
- (16) Indicator

2 Cupfiller Switch(Fig.3-5) A E

Momentarily press the cupfiller switch , water will come out from the cupfiller nozzle for 3 seconds and stops automatically.

Also the bowl flush starts and will run for 6 seconds and stops automatically. While the cupfiller is working, by momentarily pressing the cupfiller switch the cupfilling will cancel.

Also, when the cupfiller starts, the spittoon water flushes 6 seconds and stops automatically.(Synchronized Bowl Flush)

Note : The cupfiller water volume can be adjusted by the cupfiller flow control knob located inside the cuspidor body. Refer to 3-2.(3).



Fig.3-5 Cupfiller Switch

(3) Bowl Flush Switch (Fig.3-6) A E
Momentarily press the bowl flush switch, water flushes for about 5 seconds and stops automatically. (Timer mode)
Keep pressing the bowl flush over 2 sec.,water flushes continuously until the bowl flush switch is pressed again. (Continuous mode)
While the bowl flush is working by momentarily pressing the bowl flush switch the bowl flush will stop.
Note 2 : Clesta II unit can be set to timer mode (standard setting) and continuous mode for bowl flush.
The bowl flush water volume can be adjusted by the bowl flush

(4) Light Pack Switch (Optional)(Fig.3-7) Pick up the fiber optic handpiece from the holder, momentarily press the light

pack switch, and the LED illuminates in green and the fiber optic power turns on until the light pack switch is pressed again.

flow control knob located inside the cuspidor body. Refer to 3-2.(3).

(5) Chair Manual Control	Switches (Fig.3-8) A E
a. Seat Lifting	: Keep pressing 🕥 switch until the seat
	is lifted up to the desired position.
b. Seat Lowering	: Keep pressing $$ switch until the seat
	is lowered to the desired position.
c. Backrest Reclining	g : Keep pressing 🖄 switch until the backrest
	is reclined to the desired position.
d. Backrest Raising	: Keep pressing 🔿 switch until the backrest
	is raised up to the desired position.

6 Chair Auto Mode Control Switches(Fig.3-9)

a.Momentarily press the preset-1 switch (1), the chair moves to the preset 1 position and stops automatically.

Preset 2 position operated by the preset switch (2).

Note : For preset position adjustment refer to chair Manual.

b. Auto-return Operation

Fig.3-6 Bowl Flush Switch

O-LED



Fig.3-7 Light Pack Switch

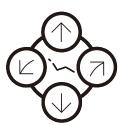


Fig.3-8 Chair Manual Control Switch

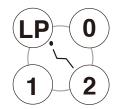


Fig.3-9 Chair Auto Mode atically Control Switch

Momentarily press (0) switch, the chair returns to the initial position (The seat is the lowest and the backrest is the upright position.) and stops automatically.

c. Last Position Memory Operation

Momentarily press (LP) switch at the reclined backrest position (treatment position), the backrest goes up to the mouth rinsing position and stops automatically. Momentarily press (LP) switch again, the backrest returns to the previous treatment position and stops automatically.

d. Emergency Stop (Safety Stop)

During automatic procedure (preset, auto-return and last position memory), momentarily pressing any chair control switch will cancel the automatic movement immediately.

(7) Safety Lock Indicator (Fig.3-10) [A][E]

The safety lock indicator illuminates amber when the safety lock device is working.

Note : Please refer to page 18 4. Safety Lock Device.

LED LOCK

Fig.3-10 Safety Lock Indicator

(8) Coolant Water ON/OFF Switch (Fig.3-11)

When a handpiece is picked up and this switch is pressed, both LED A (air) and LED W (water) lights up, the coolant water and air comes out from the handpiece. In case of air motor or air turbine, switching between spray (both of LED A and LED W are lit) and OFF occurs when this switch is pressed. In case of electric scaler, switching between water only (LED W is lit) and OFF occurs when this switch is pressed, regardless of the mode. In case of micromotor, either the 2-mode or the 4-mode can be selected by mode select setup.When this switch is pressed in the 2-mode setup, switching between spray and OFF occurs. In case of 4-mode setup, switching occurs in the sequence indicated below each time when this switch is pressed: Spray to Water only to Air only to OFF. As for the mode setting, Please refer to 3-1.(3) h.

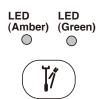


Fig.3-11 Coolant Water ON/OFF Switch



Fig.3-12 Electric Motor Speed Set Switch

1) Switching to limit rotation speed (limit mode)

The indicator indicates the selected mode.

(9) Electric Motor Speed Set Switch (Optional) (Fig.3-12)

Two different modes, limit mode and preset mode, are available as

micromotor rotation speed modes. Pressing this switch each time changes

the speed mode: Limit speed -> SET1 -> SET2 -> SET3 -> Limit speed.

Pick up the micromotor from the holder, and press this switch to select limit mode. For selecting the upper limit in the limit mode, press either plus switch or minus switch. The upper limit of the micromotor rotation speed changes in three steps (or 5 steps).

- Upper limit of rotation speed in case of 3 steps: 10000/20000/40000 min⁻¹(rpm)

- Upper limit of rotation speed in case of 5 steps: 5000/10000/20000/30000/40000 min⁻¹(rpm) See item g. on page 13 for setting in 5 steps.

The micromotor rotation speed can be varied in the range of up to the selected upper limit by sliding the foot controller pedal right or left. The rotation speed range varies by the micromotor type.

2) Switching to preset rotation speed (preset mode)

Pick up the micromotor out of the holder, press this switch, and then select preset mode (SET1 to 3). The rotation speed in this mode can be changed by pressing plus switch or minus switch. Press store switch for storing the changed rotation speed. When the foot controller is depressed upon selection of preset mode (SET1 to 3), the micromotor runs at the fixed rotation speed indicated on the indicator.

(1) Electric Motor Rotation Direction Switch (Optional) (Fig.3-13) E After picking up the electric motor from the holder, the electric


motor rotation direction can be changed by momentarily pressing this switch the rotation direction will be indicated by the amber and green LEDs.

Indicator in green : Forward Rotation

Indicator in amber : Reverse Rotation

Note: Do not change the electric motor direction while the motor is running. When the electric motor with setting Reverse Rotation is returned the holder and picked up again, a buzzer sounds.

(1) Dental Light ON/OFF Switch (Fig.3-14) Switch for on/off the dental light.

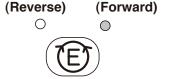

Green LED

Amber LED





Fig.3-14 Dental Light Switch

12 Function Switch (Fig.3-15)

Use this switch for setting various working conditions. As for the setup procedures, please refer to 3-1.(3) a-i.



- (3) Function Switch Setup Procedure
- a. Timer

Timer can be set maximum 90 mins. 50 secs. in 10 secs. segment.

rt Switch (Electric Motor Speed Set Switch)

1) To set a timer

Momentarily press the function switch, and set the time by pressing decrease switch and increase switch.

(+) ... Minimum setting time by switch is 10 seconds.

- Minimum setting time by switch is 1 minute.
- The setting time is indicated on the function indicator.

Momentarily press the start switch to start timer. The end of setting time is informed by electronic sounds. Example: Setting time 3 minutes 30 seconds is indicated as 03:30 in the function indicator.

2) Preset time setting

F	$\overline{}$	+		0	1	2	LP
---	---------------	---	--	---	---	---	----

Function Switch Decrease Switch Increase Switch Store Switch 0 Switch 1 Switch 2 Switch LP Switch Four preset time can be set. (0) (1) (2) (LP)

Momentarily press function switch, and set the time by pressing decrease switch and increase switch.

Press the store switch, then press the (0)(1)(2) or (LP) switch to store in memory.

3) Preset time operation

F	()	(1)	(2)	(LP)	
Function Switch	0 Switch	1 Switch	2 Switch	LP Switch	Start Switch (Electric Motor Speed Set Switch)

Press the function switch, then press the (0)(1)(2) or (LP) switch to choose desired preset number. Press start switch to start timer.

4) Cancel the timer during time countdown

$$(\mathbf{F})$$

Function Switch

Start Switch (Electric Motor Speed Set Switch)

Momentarily press function switch, then press start switch to cancel timer.

b. Group Selection Mode

The group selection mode is a convenient function for the dental clinic where a multiple (up to four groups) of dentists work with one unit. The following functions can be set for each group.

Preset position for chair. / Preset rotation speed for micro motor.

To set the group.

1) Momentarily press the function switch twice on main control panel, and the function indicator will indicate group number.

2) Momentarily press the 0,1, 2 or LP switch on main control panel to set one of 4 group.

0; Group1 / 1; Group2 / 2; Group3 / LP; Group4

c. Flush out system

The CLESTA II is equipped with two types of flush out system.

Short time flush out is for cleaning handpiece water lines.

Long time flush out is for handpiece water lines, bowl flush water line and cupfiller water line.





1) Short time flush out

Momentarily press the function switch three times and momentarily press the decrease switch. Pick up the handpieces from the holder and set them in the cuspidor bowl. By momentarily pressing the foot controller this starts short time flush out. Water comes out from the handpiece and stops automatically after 40 seconds. During flush out, by momentarily pressing any one of unit control switch or foot controller will

2) Long time flush out

cancel flush out immediately.

Momentarily press the function switch three times and momentarily press the increase switch. Pick up the handpieces from the holder and set them in the cuspidor bowl.

By momentarily pressing the foot controller this starts long time flush out for 5 minutes.

Then, cupfiller and bowl flush out starts and stops automatically in another 5 minutes.


During flush out, momentarily pressing the any one of unit control switch or foot controller will cancel flush out immediately.

* Operating instructions for flush out system (optional) |A|

Pick up the handpiece(s) from handpiece holder and set them in the cuspidor bowl.

Water flows from handpiece(s) when turn on the flush out switch.

To stop flush out operation, turn off the flush out switch during flush out.

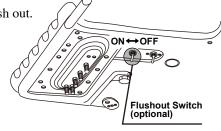

d. Control panel switching sound on/off

Pressing a switch on the control panel makes an electronic sound. This sound can be eliminated as follows;



ch Decrease Switch

+ Increase Switch


Momentarily press the function switch four times and momentarily press the decrease switch. To return to original setting.

Momentarily press the function switch four times and momentarily press the increase switch.

e. Fiber optic handpiece lighting mode (Optional)

In case that fiber optic handpiece is installed, the fiber optic turns on when the handpiece is taken out of the holder, and turns off when the handpiece is returned to the holder.

This could be changed to fiber optic turns on when the handpiece is taken out of the holder and drive air pedal of foot control is activated.

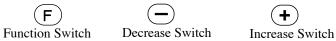



Momentarily press the function switch five times and press the decrease switch. To return to original setting.

Momentarily press the function switch five times and press the increase switch.

f. Electronic sound for timer

Electronic sound for timer can be changed.


F	0 1 2 LP
Function Switch	Chair Auto Mode Switch

Momentarily press the function switch six times.

Momentarily press one of chair auto mode switch (0,1,2,LP) then the new electronic sound is to be memorized.

g. Micro motor maximum speed setting (Optional)

The maximum rotation speed of the micro motor can be selected in 3 steps (10000,20000,40000 min⁻¹(rpm)). This function can be changed to 5 steps (5000,10000,20000,30000,40000 min⁻¹(rpm) as follows:

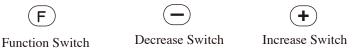

Momentarily press the function switch seven times and press the increase switch. To return to original setting.

Momentarily press the function switch seven times and press the decrease switch.

h. Coolant Water ON/OFF Switch

In case of micromotor, either the 2-mode or the 4-mode can be selected by mode select setup.

When this switch is pressed in the 2-mode setup, switching between spray and OFF occurs.

In case of 4-mode setup, switching occurs in the sequence indicated below each time when this switch is pressed: Spray to Water only to Air only to OFF.

F	$\overline{}$	+	₹¶)
Function Switch	Decrease Switch	Increase Switch	Coolant Water
			ON/OFF Switch

To set 2 mode

Momentarily press the function switch eight times and press the decrease switch.

To set 4 mode

Momentarily press the function switch eight times and press the increase switch.

i. Cupfiller and bowl flush

Cupfiller and bowl flush are set to operate together (when the cupfiller switch is activated, bowl flush also starts).

To make these operate independently.





Increase Switch

Momentarily press the function switch nine times and press the decrease switch. To return to original setting.

Momentarily press the function switch nine times and press the increase switch.

(4) Doctor Table Chassis (Fig.3-16)

A.The handpiece spray water control knobs are located under the doctor table.

Each handpiece spray water control knob is marked 1-4 from the left side HP1,HP2,HP3,...

The handpiece spray water volume can be controlled independently.

- *The HP4 is optional.
- B. Doctor's Syringe Flow Control Knobs (Fig.3-16)

Doctor's syringe flow control knobs are located under the doctor table.

The flow control knobs adjust the doctor's syringe air and water flow volume.

The yellow capped knob is the air flow control knob,

the blue capped knob is the water flow control knob.

- Note : Turning the control knob counterclockwise will increase the flow volume and turning clockwise will decrease.
- C. Handpiece Air Pressure Gauge (Fig.3-16)During a handpiece is working, the handpiece drive air pressure is shown in the pressure gauge.
- (5) Balance Arm Air Brake (Fig.3-17)

When the master switch is ON, the balance arm is locked by the air brake.

Hand the handle and press the air brake release button to adjust the table height. At the desired table position, release the air brake release button and the balance arm is locked.


(6) Dental Size Film Viewer (Fig.3-18)

Film viewer ON/OFF switch is located on right-hand side of the film viewer.


Press the switch and the film viewer turns on. Press again and the film viewer turns off.

Handpieces

Refer to handpiece manufacturers operating instructions.

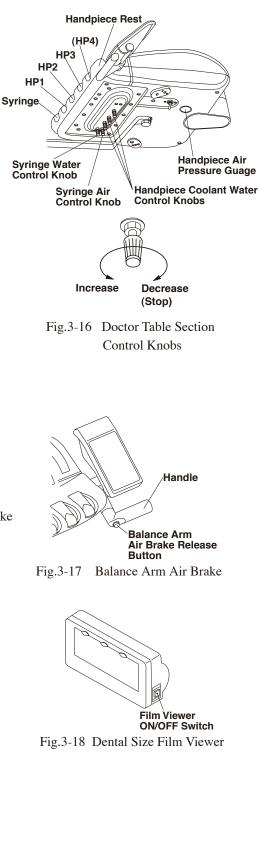

3-2. CUSPIDOR UNIT SECTION

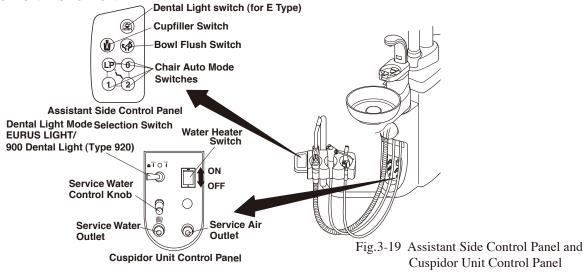

(1) Assistant Side Control Panel (Fig.3-19)

Cupfiller switch, bowl flush switch, dental light switch and chair auto mode switches are located on the assistant side control panel. Refer to 3-1.(2) Main Control Panel

- (2) Cuspidor Unit Control Panel (Fig.3-19)
 - A. Dental Light Mode Selection Switch
 - 900 Dental Light (Type 920) / EURUS LIGHT

Dental light can be operated (ON/OFF) either by the touchless sensor switch located

on the light head or by the manual switch on the cuspidor unit control panel.

To operate by the touchless sensor switch --- Set the switch lever to SENSOR Side

To operate manually ----- Set the switch lever to Manual Side

Set the switch lever to Centre for OFF

B. Water Heater Switch (Fig.3-19)

Water heater switch is located on the cuspidor unit control panel. Turn on the water heater switch and the cupfiller water will warm up.

C. Service Water Outlet (Fig.3-19)

The service water outlet provides a quick-connection for water.

D. Service Water Outlet Control Knob (Fig.3-19)

The water volume from the service water outlet can be adjusted by the service water outlet control knob. Turning the knob counterclockwise will increase the flow volume and turning clockwise will decrease.

E. Service Air Outlet (Optional) (Fig.3-19)Service air outlet provides a quick-connection for air.

(3) Cuspidor Unit Body

A. Assistant's Syringe Flow Control Knobs (Fig.3-20) Assistant's syringe flow control knobs are located in the cuspidor unit body.

The yellow capped knob is to adjust the assistant's syringe air flow volume, and the blue capped knob is to adjust water flow volume.

B. Cupfiller Flow Control Knob (Fig.3-20)
Cupfiller flow volume can be controlled by the cupfiller flow control knob. (Pinch valve system)
Loosen the lock nut and adjust cupfiller water flow volume by turning the knob.
Tighten the lock nut after adjustment.

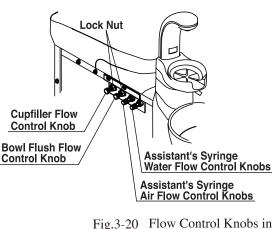


Fig.3-20 Flow Control Knobs in Cuspidor Unit Body C. Bowl Flush Flow Control Knob (Fig.3-20)


Bowl flush flow volume can be controlled by the bowl flush flow control knob. (Pinch valve system) Loosen the lock nut and adjust bowl flush water flow volume by turning the knob.

Tighten the lock nut after adjustment.

- Note : Turning a knob counterclockwise increase flow volume and turning clockwise will decrease.
- (4) Assistant Instrument Holder (Fig.3-21)

When picked up an instrument (Saliva ejector or Vacuum handpiece) from the assistant holder this starts the instrument working automatically.

Returning the instrument to the holder stops automatically. Saliva ejector handpiece and vacuum handpiece have stop valves to close and adjust suction power.

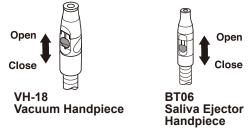




(5) Sensor Cupfiller (Optional) (Fig.3-22)


Place the cup (paper cup) on the cupfiller base and water comes out from the cupfiller nozzle fills up the cup and stops automatically.

When cupfiller starts, the bowl flush also starts, and will run for about 6 sec. and it stops automatically. While filling the cup, momentarily press the cupfiller switch, and this will cancel the cupfiller water flow. During bowl flush, momentarily press the bowl flush switch and this will cancel the bowl flush water flow. Note: Use only suitable disposable paper cup (dental paper cup).

Use only an empty cup. A cup with some water in could cause an over flow.

The sensor cupfiller needs over 2 seconds interval between cup filling.

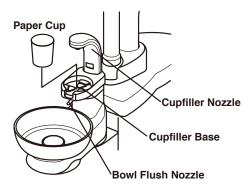


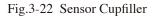

(6) Cuspidor Bowl Rotation (Optional) (Fig.3-23)

The cuspidor bowl can be rotated 90° each. (inside & outside)

Do not hold the cupfiller nozzle when rotating the cuspidor bowl; be sure to hold the bottom of the spittoon bowl. Rotating via the cupfiller nozzle may result in damage to the cupfiller nozzle.

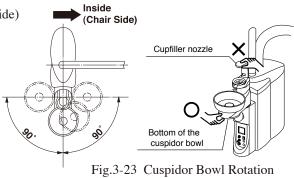

- (7) Clean Water System (Optional) (Fig.3-24, Fig.3-25)
 - The water bottle CITY/BOTTLE select switch can be changed between municipal water and water bottle. CITY.... Tap water
 - **BOTTLE . . . Water bottle**
 - When the pressurization ON-OFF switch is set to ON, the water bottle may be used. (Fig.3-24)
 - The water bottle pressure control knob adjusts the pressure of water bottle. The pressure increase when the knob is turned

clockwise, and decreases when the knob is turned counterclockwise. (Fig.3-24)

Remove or attach the bottle (Fig.3-25) Turn off the pressurization ON-OFF switch. Water bottle may be remove by turned in direction A.

It is attached when turned in direction B.

- The water tank is intended only for use with purified water, distilled water and pure water.
- Do not use mouthwash or electrolyzed water, as they may cause clogged tubing or affect internal valves and equipment.
- Adjust the air supply pressure for the water tank to 200 kPa or less. An excessively high pressure may cause damage to the water tank.


(8) Height Adjustable Assistant Arm (Optional) (Fig.3-26) Press the lock release button and

raise arm to adjust the assistant

holder height.

- Position at desired height and
- release the lock button after that.
- Note: Support the arm with your hand until it is positioned at the desired height.

Water Bottle Pressure Control Knob CITY/BOTTLE Select Switch Water Bottle Water Bottle Pressure Control Knob

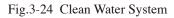

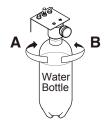


Fig.3-25 Remove or Attach the bottle

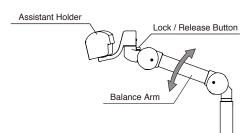


Fig.3-26 Height Adjustable Assistant Arm

(9) Dental Light (Fig.3-27)

Please refer to operating instruction for dental light.



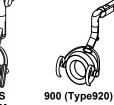






Fig.3-27 Dental Light

4-1 T :-1-4

(10) Monitor Bracket (Optional) (Fig.3-28) PC monitor can be mounted on this bracket.

Note: Do not apply an excessive load or shock to the monitor or monitor bracket.



Fig.3-28 Monitor bracket

3-3. FOOT CONTROL SECTION

- (1) Foot Control (Type A2) (Fig.3-29 & 3-30) A
 - A. Drive Air Pedal

Pick up a handpiece from the handpiece rest (pull the rod about 10 degrees to forward) and depress the drive air pedal, the handpiece starts running.

- B. Coolant Water Switch Coolant water switch allows handpiece coolant water to be turned on or off.
- C. Chip Blower Pedal

By depressing the chip blower button, the chip air will come out from handpiece without the bur rotating.

- (2) Electric Motor Foot Control (Type SE) (Fig.3-31)
 - A. Drive Air Pedal

Pick up a handpiece from the handpiece rest (pull the rod about 10 degrees to forward) and depress the drive air pedal, the handpiece starts running.

B. Coolant Water Switch

Momentarily depressing the coolant water switch is changed handpiece coolant water and air situation.


The situation is shown on the main control panel.


Refer to page 10, (8) Coolant Water ON/OFF Switch.

C. Chip Blower Button

By depressing the chip blower button, the chip air will come out from handpiece without the bur rotating.

D. Electric Motor Rotation Control (Optional)

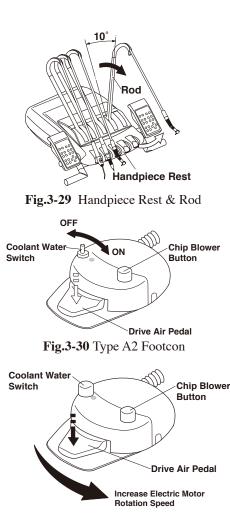

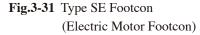

Pick up the electric motor from the handpiece holder and while pressing down slide drive air pedal horizontally to right, and the electric motor will start running.

The rotation speed increases by sliding the drive air pedal further to the right.

The speed control by the foot control is within the limits of the electric motor speed setting.

E. Coolant Water Switch / Electric Motor Rotation Direction Switch (Optional)

The coolant water switch can be changed for electric motor rotation direction switch.

To change to electric motor rotation direction switch.

Keep depressing the coolant water switch until buzzer sound (about 2 sec.).

To return to original (coolant water switch) setting.

Keep depressing the coolant water switch until buzzer sound (about 2 sec.).

The electric motor rotation direction is indicated on the control panel by LED.

Please see page 10, 10 Electric Motor Rotation Direction Control Switch

- (3) Foot Switch (Optional) (Fig.3-32)
 - A. Chair Manual Control Switches
 - B. Chair Auto Mode Switches

Note: Please refer to page 9 (5) (**Fig3-8**)(**Fig.3-9**).

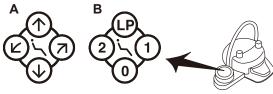


Fig.3-32 Foot Switch

4. SAFETY LOCK DEVICE

In the following cases the safety lock device to lock the chair movement is activated.

- 1. When the pedal of the foot controller is depressed.
- 2. Press any chair operation switches for stopping the automatic movement.
- 3. During setting with the function switch on the doctor control panel.

4. When the cuspidor bowl is rotated toward the patient side. (Pedestal type)

Note: Please refer to page 9, (7) (Fig3-10).

5. CARE AND MAINTENANCE

Cleaning Unit

Cleaning and Sterilization of product exterior

Clean the metallic parts with a dry soft cloth. Wipe off water immediately if water is put on the product. Water may cause rusting.

Clean the resin parts with a wet soft cloth.

Use FD333 or FD366 made by Durr or ethanol for cleaning and bacteria elimination from the product exterior.

- When the surface of the control panel is cleaned with disinfectant, etc., wipe off disinfectant completely. If it penetrates into the back of the sheet, the membrane switches may malfunction.
 Wipe off water and residual disinfectant immediately. This could cause corrosion, damage or
- incorrect operation of the unit.

• Immediately wipe off any water spills or leakage on the floor. This could cause damage to the product, decreased strength of the floor may lead to physical injury including fall, or property damage.

(1) Drain Cap, Basket Strainer, Cupfiller base, and Cuspidor Bowl

- A. Detach the drain cap, and clean the basket strainer.
 Remove the cupfiller base and clean it.
 Turn the bowl flush nozzle counterclockwise, avoiding the clash with the cuspidor bowl. (Fig.5-1)
- B. Turn the cuspidor bowl to left and right and remove it upward. (Fig.5-2)

- Turn off the main switch before removing the cuspidor bowl.
- Cuspidor bowl is made of glass or ceramic and may be broken when given the impact. When and after removing or attaching the cuspidor bowl, beaware of the handling, avoiding the impact, hitting and falling.

When carrying the cuspidor bowl, hold it firmly with both hands while placing one hand at the bottom of the cuspidor bowl.

- Do not clean the cuspidor bowl using hot water because it may be broken.
- Put on the durable gloves when cleaning.
- C. When it is hard to attach the cuspidor bowl after cleaning, wet the inserting surface for an easy attachment.After the attachment, make sure the cuspidor bowl is attached firmly. (Fig.5-3)
- D. After attaching the cuspidor bowl, return the bowl flush nozzle to the original position and attach the cupfiller base.

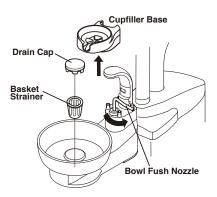


Fig.5-1 Drain Cap, Basket Strainer, and Cupfiller Base

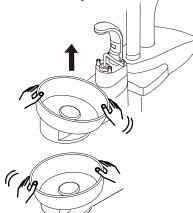


Fig.5-2 Remove the Cuspidor Bowl

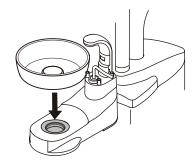


Fig.5-3 Attach the Cuspidor Bowl

(2) Solid Collector (Fig.5-4)

Pull out the solid collector filter and clean it. After cleaning, refit the solid collector firmly.

(3) Tubings and hoses

Tubings and hoses can be cleaned with FD333 or FD366 made by Durr.

(4) Cleaning vacuum and saliva ejector lines (Fig.5-5)

The sucking unit comes into contact with secretions, spit and blood that contain bacteria every day. Be sure to clean and sterilize it at the end of each work day. Recommended cleaner : Orotol Plus made by DURR.

Do not use strongly acidic cleaning agents or alkaline pipe cleaning agents, which may cause corrosion of metals, etc/

(5) Vacuum hose and saliva ejector hose are detachable from the cuspidor unit (**Fig.5-6**)

Turn the hose connector 90° counterclockwise. The hose can be removed from the cuspidor unit.

Insert the hose connector and turn through 90° clockwise to reconnect.

Note: After daily operation, run two cups of clean water through handpieces to clean inside.



Fig.5-4 Solid Collector Filter



Fig.5-5 Cleaning Vacuum/ Saliva Ejector Lines

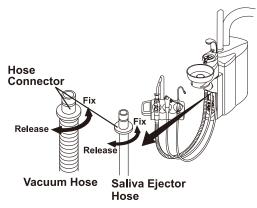


Fig.5-6 Vacuum and Saliva Ejector Hose

- (6) Vacuum Handpiece and saliva ejector Handpiece
 - * Cleaning amd sterilization must be done every after use to patients.

For effective sterilization, washing for removing contamination and immersion by a cleaning agent are required. Then, rinse by water in order to remove residual procedures from cleaning to sterilization.

* Use the disposable saliva ejector tip.

Disassembly (Fig.5-7)

Disassemble the handpiece for the preparation of cleaning as the following figures show. Pull the hose connector to disconnect the vacuum hose.

Cleaning by hand

- A. Wipe off the surface contamination by a cloth while rinsing the surface by running clean warm at 40±5 degrees. (Fig.5-8) Scrub the intubation or hole, slide groove and filter by a cleaning brush or by a tooth brush with running clean warm water at 40± degress. (Fig.5-9) Wipe off by a cloth for the area which brush is unable to reach. (Fig.5-10)
- B. Check whether contamination is removed or not after cleaning.

Continue the cleaning if contamination is remained.

- C. Immersed with an alkaline disinfection or detergent for 5 minutes. (We recommend to use ID212 made by DURR)
- D. Rinse thoroughly by distilled water at ordinary temperature or by clean water for more than 1 minute.

Cleaning must be done within 1 hour after use. Replace with a new vacuum handpiece and saliva ejector handpiece for following cases.

- Any waste material can not be removed by clogged hole.
- Contamination and solid material attached to vacuum handpiece and saliva ejector can not be removed.

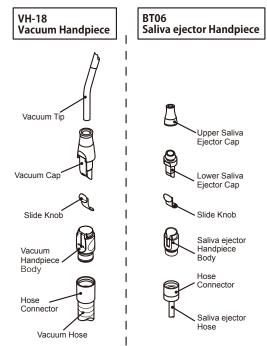


Fig.5-7 Vacuum Handpiece and Saliva Ejector Handpiece



Fig.5-8 Surface cleaning

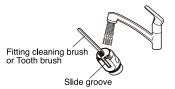


Fig.5-9 Cleaning of hole and sliding part



Fig.5-10 Cleaning the area (brush is unable to reach)

Assembly before sterilization (Fig.5-11)

Vacuum Tip/Vacuum Handpiece/Saliva Ejector Handpiece can be autoclave. Vacuum handpiece and saliva ejector handpiece have to assemble before autoclave.

When assembling the vacuum handpiece, align the two surfaces; a flat surface of the cap and the surface of the sliding groove of the body. Then, slide the cap into the body straight.

* Do not insert it twisted.

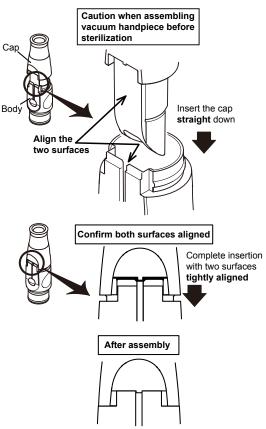


Fig.5-11 Assembly before sterilization (Vacuum Handpiece)

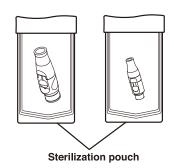


Fig.5-12 Vacuum Handpiece and Saliva Ejector Handpiece

Sterilization (Fig.5-12)

A. Insert the handpiece in a sterilization pouch and seal it.B. Autoclave for 3 min. at 134°C and dry for 15 min.Sterilization with autoclave is pemitted up to 250 times.



- Sterilization by class B cycles.
- Sterilization temperature is 135°C or less.
- The cap and body are made of resin. They may become deteriorated if they are sterilized in an autoclave many times.
- After autoclave sterilization, the cap, body, and slide knob are subject to discoloration, which does not have a negative effect on performance.
- The slide knob can be autoclave100 times and is expendable supplies.
- Dry naturally if the temperature for drying process is to exceed 135°C.
- If damage occurs to the Sterilization pouch, discard and sterilize again using a new pouch.

Storage

After sterilizing the handpiece, keep it in the sterilization pouch and store in a dark and cool place.

(7) Belmont 77 Syringe

* Cleaning amd sterilization must be done every after use to patients.

Disassembly (Fig.5-13)

Remove the nozzle from syringe by turning it in direction A.

Cleaning by hand

- A. Wipe off the surface contamination by a cloth while rinsing the surface by running clean warm water at 40±5 degrees. Scrub the tip and joint part of nozzle by a cleaning brush or by a tooth brush with running clean warm water at 40±5 degrees. (Fig.5-14)
- B. Check whether contamination is removed or not after cleaning. Continue the cleaning if contamination is remained.
- C. Immersed with an alkaline disinfection or detergent for 5 minutes. (We recommend to use ID212 made by DURR)
- D. Rinse thoroughly by distilled water at ordinary temperature or by clean water for more than 1 minute.

Cleaning must be done within 1 hour after use. Replace with a new syringe nozzle if contamination and solid material attached to syringe nozzle can not be removed.

Sterilization (Fig.5-15)

The nozzle can be sterilized with autoclave.


A. Insert the handpiece in a sterilization pouch and seal it.


B. Autoclave for 3 min. at 134°C and dry for 15 min.

Sterilization with autoclave is permitted up to 250 times.

Apply Vaseline thinly and evenly to the two O-rings after sterilization.

Tighten the nut firmly in direction B to undo the nozzle.

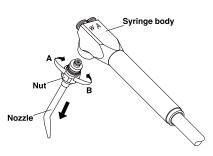

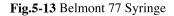

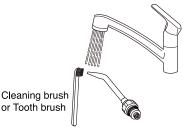


Fig.5-14 Cleaning Belmont 77 Syringe Nozzle



Sterilization pouch

Fig.5-15 Sterilization Belmont 77 Syringe Nozzle

- Sterilization by class B cycles.
- Sterilization temperature is 135°C or less.
- Dry naturally if the temperature for drying process is to exceed 135°C.
- If damage occurs to the Sterilization pouch, discard and sterilize again using a new pouch.
- Before use, make sure that the nut is firmly tightened.

Storage

After sterilizing the nozzle, keep it in the sterilization pouch and store in a dark and cool place.


(8) Handpieces

Observe the instructions given in the package insert and Instruction Manual included with the handpiece to clean it.

(9) Air Filter Drain Valve (Fig.5-16)

Air filter drain valve is located under the maintenance lid.

Once a week open the drain valve and drain off water from the air line.

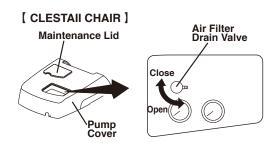

[CLESTAII CHAIR (EURUS TYPE)]


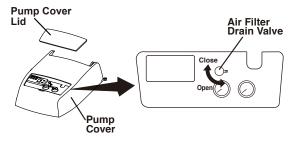


Fig.5-16 Air Filter Drain Valve

[CLESTAII CHAIR]

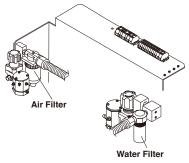

[CLESTAII CHAIR (EURUS TYPE)]

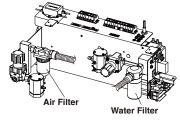


Fig.5-17 Filter

(10) Filter Replacement (Fig.5-17)

The water filter in the junction box needs to be replaced at least once a year. The air filter in the junction box needs to be

The an inter in the junction box needs to be

replaced at least once every three years.

Contact your local service representative for replacement.

6. MAINTANANCE AND INSPECTION

6-1. Guide for daily maintenance and inspection (Maintenance and inspection by user)

- •Management of maintenance and inspection of medical equipment should be implemented by the user (medical institution).
- For safe use of this product, it is necessary that inspection should be conducted in the specified frequency on the tems described below.

No.	Item	Frequency	Inspection method and diagnosis	Influence if inspection not conducted	Maintenance required in case of nonconformity
1	Check of safety functions	Before start	 Make sure the chair movement stops by any of the following actions. (1) when foot controller pedal is depressed. (2) During chair auto movement , press for any chair operation switch. (3) While setup is in progress with function switch on the doctor membrane switch panel. (4) When the spittoon bowl is turned to patient side (Pedestal type) 	Unexpected personal injury and troubles may arise due to motion of the chair during medical treatment and due to pinching between doctor section and chair.	Contact your dealer or our office if any abnormality arises.
2	Check for leakage of water and air	Before start	Leakage of water and air shall not be observed around the product.	The product will not normally work, and troubles may arise.	Contact your dealer or our office if any abnormality arises.
3	Cupfiller	Before start	When a paper cup is placed on the cupfiller, the cup shall be detected and cupfilling shall be executed. * Malfunction may arise if the cup is of another material grade (such as stainless steel and plastics) or if the paper cup is of dark color or pattern.	Cupfilling may not be executed.	Execute re-inspection in accor- dance with "Method for operation" described in the instruction manual. Contact your dealer or our office if recovery is not achieved as a result of re-inspection.
4	Check of motions of equipment	Before start	 Air turbine revolution, water flow, air flow and so forth shall be free of abnormality. Micromotor revolution, water flow and so forth shall be free of abnormality. Scaler vibration, water flow and so forth shall be free of abnormality. 	Troubles such as injury in patient's oral cavity and equipment failure may arise.	Control the water flow in accor- dance with "Control of compo- nents" described in the instruction manual. If any other abnormality arises, refer to the instruction manual attached to individual equipment. Contact your dealer or our office if recovery is not achieved.
5	Check of air turbine bar	For each patient	Appropriate bar shall be positively mounted. Make sure to refer to the instruction manual attached to individual equipment.	The bar will not normally work and troubles may arise.	If abnormality such as flaw and deformation is found on the bar, replace the bar in accordance with the instruction manual attached to individual equipment.
6	Check of scaler tip	For each patient	Appropriate tip shall be positively mounted and be correctly used. Make sure to refer to the instruction manual attached to the scaler.	The tip will not normally work and troubles may arise.	If the tip was worn or deformed, replace the tip in accordance with the instruction manual attached to the scaler. Contact your dealer or our office if any other trouble arises.
7	Check of tightness of syringe nut	For each patient	The nut for fixing the nut of Type 77, 3-way syringe shall be positively tightened.	Troubles may arise if the nut comes off.	Turn and positively retighten the nut that fixes the nozzle.
8	Matters attached to micromotor	After closing	Excessive handpiece oil or the like shall not be attached to the motor section.	The motor section will not work normally and troubles may arise.	Execute care in accordance with the instruction manual attached to individual micromotor.

No.	Item	Frequency	Inspection method and diagnosis	Influence if inspection not conducted	Maintenance required in case of nonconformity
	Care Vacuum and saliva ejector handpiece	After closing	Flush the suction line, and then clean the filter of the vacuum or saliva ejector handpiece.	Faulty suction may arise.	Clean the suction line and filter in accordance with "Method for care" described in the instruction manual.
10	Care Cuspidor section	After closing	Clean the cuspidor and dust filter.	Faulty water drainage may arise.	Execute dust removal and cleaning in accordance with "Method for care" described in the instruction manual.
11	Care Solid collector	After closing	Clean the filter of the solid collector.	Vacuum suction will become weak.	Clean the filter in accordance with "Method for care" described in the instruction manual.
12	Care Exterior	After closing	Chemical, filthy water and so forth shall not be found (attached or remaining) on the product exterior.	Discoloration and deterioration to the exterior, and corrosion and rusting to metallic components may arise.	Execute wiping in accordance with "Method for care" described in the instruction manual.
13	Check of main switch and main valves	After closing	The product main switch shall be off, and air main valves shall be closed.	Product failure and troubles may arise.	Contact your dealer or our office if the main switch cannot be turned off or if the main valve cannot be closed.
14	Product's moving parts	Once every week	No abnormal noise or the like shall be produced from product's moving parts when the product is operated.	The product will not normally work and troubles may arise.	Contact your dealer or our office if any abnormality arises.
15	Care Drain valve	Once every week	Water may enter the air line, and equipment failure may arise.	Drain the water from the air filter drain valve.	Drain the water from the air filter in accordance with "Method for care" described in the instruction manual.
16	Check of water pressure and air pressure	Once every month	Check the water pressure and air pressure by reading pressure gauges in the U-Box section. Basic set pressure : Water: 0.1 to 0.2 MPa Air: 0.45 to 0.5 MPa	The product will not normally work, and troubles may arise.	Contact your dealer or our office if the set pressure is abnormally high or low.
17	Check of conditions of table section	Once every month	The table shall be free of inclination, and water shall not flow on the table.	Injury caused by falling of goods located on the table and other troubles may arise.	Contact your dealer or our office if any abnormality arises.
18	Care Oil mist separator	Once every month	The oil level in the oil mist separator shall be lower than the red line.	Normal output will not be produced due to inferior handpiece exhaust.	Discharge the oil in accordance with "Method for care" described in the instruction manual.

6-2. Guide for Periodical Check-up

- Some parts and components of the products are degraded or deteriorated depending on the frequency of use. Annual check-up and maintenance, as well as replacement of consumable parts, are required.
- The required parts (including consumable parts) are listed below. It may be different from the following list depending on the option of the unit.
- For check-up and repair, call a technician of our authorized dealer.

Parts and components that require periodical check-up

No.	Parts Description	Standard Lifetime	No.	Parts Description	Standard Lifetime
1	Vacuum handpiece body	3 years	8	Regulator	3 years
2	Saliva ejector handpiece body	3 years	9	Valves	3 years
3	Foot controller	5 years	10	Switches	5 years
4	Water supply hose	3 years	11	Film viewer body part	5 years
5	Drain hose	3 years	12	Pressure gauge	3 years
6	Air supply hose	3 years	13	Arm section of moving part	7 years
7	Electric wiring of moving part	5 years	14	Control PCBs.	5 years

Consumable parts

No.	Parts Description	No.	Parts Description
1	Valve for vacum handpiece body	6	Filter for oil mist separator
2	Vacuum tip	7	Filter (Air & Water)
3	Handpiece tubings	8	O-ring, Packing, Diaphragm
4	Vacuum hose		
5	Saliva ejector hose		

AWARNING

Execute the maintenance in accordance with this instraction manual and operating manual attached to each individual equipment (Dental light, Handpiece, etc..). Failure to maintain this product may lead to physical injury or property damage.

7. Detachable Parts

PARTS	PARTS
Handpiece	Saliva ejector Handpiece
Handpiece hose	Vacuum tip
Cuspidor bowl	Solid collector lid
Drain cap	Solid collector filter
Basket strainer	Male connector for water service outlet
Cupfiller base	Male connector for air service outlet
Cupfiller nozzle	Oil mist separator
Vacuum Handpiece	

8. ELECTROMAGNETIC COMPATIBILITY(EMC)

This product complies with EMC Standard EN60601-1-2:2015.

1. Precautions regarding EMC and compliance with accompanying documents

Medical electrical equipment requires special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in this booklet.

2. Effects of RF communication devices

Portable and mobile RF communication devices can affect medical electrical equipment.

3. Installation exclusion environment

Hospitals except for near active HF SURGICAL EQUIPMENT and the RF shielded room of an

ME SYSTEM for magnetic resonance imaging, where the intensity of EM DISTURBANCES is high.

4. Electromagnetic emission declaration

Guidance and manufacturer's declaration – electromagnetic emissions					
The CLESTA II Unit is intended for use in the electromagnetic environment specified below. The customer or user of the CLESTA II Unit should ensure that it is used in such an environment.					
Emissions test					
RF emissions CISPR 11	Group 1	The CLESTA II Unit only uses RF energy for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference with nearby electronic equipment.			
RF emissions CISPR 11	Class B	The CLESTA II Unit is suitable for use in all establishment including domestic establishments and those directly conne			
Harmonic emissions IEC 61000-3-2	Class A	to the public low-voltage power supply network that supplies buildings used for domestic purposes.			
Voltage fluctuations/ Flicker emissions IEC 61000-3-3	Complies				

Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

5. Electromagnetic immunity declaration 1

Guidance and manufacturer's declaration – electromagnetic immunity

The CLESTA II Unit is intended for use in the electromagnetic environment specified below. The customer or user of the CLESTA II Unit should ensure that it is used in such an environment.

user of the CLESTAIL Unit should ensure that it is used in such an environment.					
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance		
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Flooring should be wood, concrete, or ceramic tiles. If the floor is covered with synthetic material, the relative humidity should be at least 30%.		
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	The mains power quality should be that of a typical commercial or hospital environment.		
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	The mains power quality should be that of a typical commercial or hospital environment.		
Voltage dips, short interruptions and	0% UT ; 0.5 cycles	0% UT ; 0.5 cycles	The mains power quality should be that of a typical commercial or hospital		
voltage variations on power supply	0°,45°,90°,135°,180°, 225°, 270° and 315°	0°,45°,90°,135°,180°, 225°, 270° and 315°	environment. If the user of the CLESTA II Unit requires continued		
input lines IEC 61000-4-11	0% UT ; 1cycle and	0% UT ; 1cycle and	operation during mains power interruptions, it is recommended that the		
	70% UT ; 25/30 cycles at 0°, single phase	70% UT ; 25/30 cycles at 0°, single phase	CLESTAII Unit be powered from an uninterruptible power supply or a battery.		
	0%UT ; 250/300 cycles	0%UT ; 250/300 cycles			
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.		
Note U_T is the AC mains voltage prior to the application of the test level.					

6. Electromagnetic immunity declaration 2

Guidance and manufacturer's declaration – electromagnetic immunity					
The CLESTA II Unit is intended for use in the electromagnetic environment specified below. The customer or user of the CLESTA II Unit should ensure that it is used in such an environment.					
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance		
Conducted RF IEC 61000-4-6	3V 0.15MHz~80MHz 6V 0.15MHz~80MHz in ISM and amateur radio bands	3V 0.15MHz~80MHz 6V 0.15MHz~80MHz in ISM and amateur radio bands	Warning: Portable RF communications equipment (including peripherals such as antenna cablesand external antennas) should be used no closer than 30 cm(12 inches) to any part of the CLESTA II, including		
Radiated RF IEC 61000-4-3	3V/m 80MHz~2.7GHz 80% AM (1 kHz)	3V/m 80MHz~2.7GHz 80% AM (1 kHz)	cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.		
Near electromagnetic field caused by RF wireless communication devices IEC 61000-4-3	See the table below	See the table below			

7. Essential performance

Unless operated by the chair control switch, the chair will not make any movements. Unless operated by the foot controller, the handpiece will not move except for sounding a buzzer and switching the indicator on/off. Loss or decline of essential performance may cause the chair or handpiece to move unexpectedly, causing harm to the patient, operator or people or objects around the patient or operator.

Portable RF communications equipment (including peripherals such as antenna cablesand external antennas) should be used no closer than 30 cm(12 inches) to any part of the CLESTA II, including cables specified by the manufacturer.

Otherwise, degradation of the performance of this equipment could result.

Test frequency (MHz)	Modulation	IEC 60601 test level	IEC 60601 compliance level	
385	Pulse modulation ^{a)} 18Hz	27V/m	27V/m	
450	Frequency modulation ±5kHz shift 1kHz sine wave	28V/m	28V/m	
710 745 780	Pulse modulation ^{a)} 217Hz	9V/m	9V/m	
810 870 930	Pulse modulation ^{a)} 18Hz	28V/m	28V/m	
1720 1845 1970	Pulse modulation ^{a)} 217Hz	28V/m	28V/m	
2450	Pulse modulation ^{a)} 217Hz	28V/m	28V/m	
5240 5500 5785	Pulse modulation ^{a)} 217Hz	9V/m	9V/m	
Note a) The	Note a) The carrier is modulated by a square wave with a 50% duty cycle.			

Near electromagnetic field caused by RF wireless communication devices

9. HANDPIECES THAT CAN BE USED WITH THIS PRODUCT

	DESCRIPTION
~ .	LUZZANI (3-way) Minilight w/Light
Syringe	LUZZANI (6-way) Minilight w/Light
	DCI (3-way) SYRINGE (3459 Quick clean continental syringe)
Turbine	NSK Ti-Max X
Air motor	NSK EX-203 / EX-6
Micromotor	BIEN AIR MX2 / DMX2 PRO PCB. / PM 1:1
	BIEN AIR MCX / DMCX PCB. / PM 1:1
	NSK NLX plus
	NSK NLX nano
	NSK NBX / iMD PCB. / EX-6 (Only E-Type Holder / Continental Type)
Scaler	SATELEC NEWTRON LED / SP4055 NEWTRON Module + LED Drive Board
Sealer	NSK VARIOS 170 w/light
	XINETIC SCALER
	EMS NO PAIN SCALER

IMPORTANT

The connectors of our handpiece hose are designed and manufactured in accordance with ISO9168:2009 (DENTISTRY-HOSE CONNECTORS FOR AIR DRIVEN DENTAL HANDPIECES).

However, there still could be a case that the connectors of turbines or air motors may not fit into some handpieces due to the manufacturing tolerances.

Have your local authorized Belmont dealer check the connectability before purchasing the handpiece.

Except for our recommended handpieces, we shall not be liable for any problems deriving from bad connectability or their performance.

10. CHAIRS THAT CAN BE USED WITH THIS PRODUCT

- CLESTA [] CHAIR
- CLESTA II CHAIR (EURUS TYPE)
- CLAIR CHAIR

11. DENTAL LIGHTS THAT CAN BE USED WITH THIS PRODUCT

- EURUS LIGHT
- 900 Dental Light (Type 920)
- 300 Dental Light (Type 320S)
- AL-720S / AL-720M

12. ACCESSORIES

It may be different from the following list depending on the specification of the unit.

Accessories	intended use, etc.
Cotton containers	The container to take in cotton and swab for dental treatment.
Waste receptacle	Receptacle to do away with used cotton,etc.
Foot control handle	Grip for portability of foot control to reduce hand contact.
Quick connector for water service outlet Model number: MCL-04NH-1B	Quick connector that provides external equipment with water.
Quick connector for air service outlet ໜື່ມ Model number: MC-04PH	Quick connector that provides external equipment with air.
Flushout tray	The flushout tray holds the handpieces and allows the hands-free flushing.

Declaration of Conformity

We hereby declare that the product listed below complies with the essential requirements of the Medical Device Directive: 93/42/EEC and RoHS Directive: 2011/65/EU based on category 8 of Annex I.

Product : DENTAL UNIT (CLASS IIa)

Model : CLESTAII

" CLASS IIa " has been defined by the rule 9 of MDD Annex IX.

The product has been designed and manufactured in accordance with the European standards as the listed in the Declaration of Conformity.

Our full quality assurance system has been approved to MDD ANNEX II, excluding section 4 (Registration No. HD 60143723 0001) by Notified Body, TUV Rheinland LGA Products GmbH (NO.0197).

NOTE

EC REP

TAKARA COMPANY EUROPE GmbH

Berner Strasse 18, 60437 Frankfurt am Main, Germany Tel:+49 69 506878 0 Fax:+49 69 506878 20



TAKARA BELMONT CORPORATION2-1-1,Higashishinsaibashi ,Chuo-ku,Osaka, 542-0083 JapanTEL. : +81 6 6213-5945FAX : +81 6 6212 3680

BOOK NO. FEFA27Y0 Printed in Japan 2019-12